
waterpik®

Waterpik® Sonic-Fusion® 2.0 Models SF-03 | SF-04

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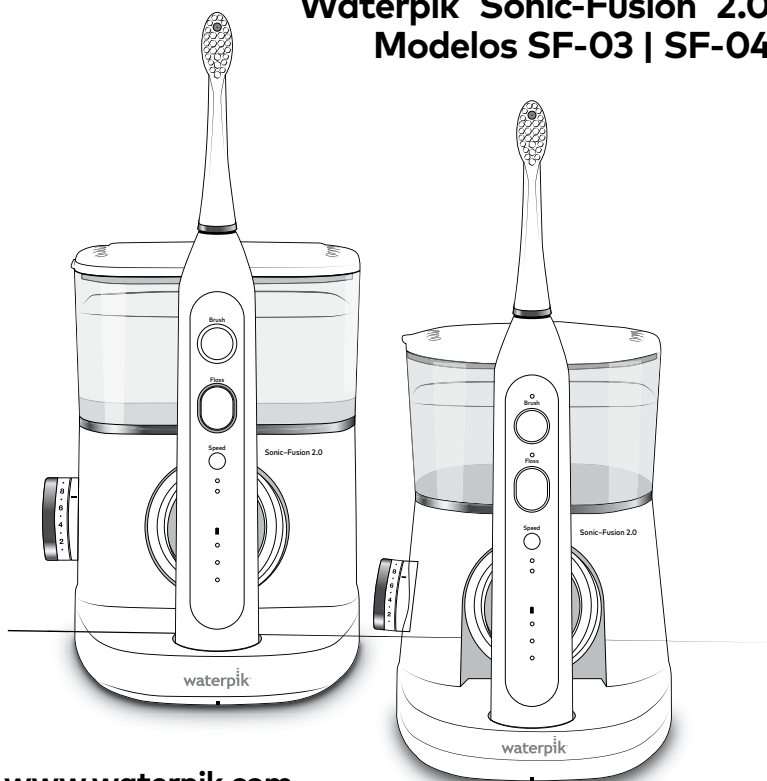IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.


- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Do not handle plug with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.

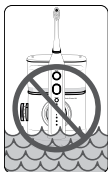WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use flossing brush heads and accessories recommended by Water Pik, Inc.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Fill reservoir with water or other dental professional recommended solutions only.
- Remove any oral jewelry prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.

- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist and/or medical professional before using this product or any other oral hygiene aid.
- Do not use iodine or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet. Contact your doctor with any questions.
- Do not use the toothbrush motor for more than 5 minutes in each 10 minute period.
- This product does not contain any user serviceable parts. Do not attempt to open.
- Using toothpastes containing peroxide may damage the plastic housing of this product.
- Using abrasive or chemical cleaners may damage the plastic housing of this product.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- The supply cord cannot be replaced. If the cord is damaged the appliance should be scrapped.
- Do not use this product if it has a damaged cord or power supply, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Water Pik, Inc. or visit www.waterpik.com.
- **IMPORTANT** This device uses wireless communication to control the ON/OFF function of the water flosser pump motor from the toothbrush power handle. To avoid electromagnetic interference, do not operate this device in restricted locations such as aircraft or specially marked areas in hospitals.


SAVE THESE INSTRUCTIONS

TABLE OF CONTENTS

Product Description	5
Preparing The Unit	6
Brushing and Flossing	6
Care and Maintenance	8
Limited Three-Year Warranty	11

PRODUCT DESCRIPTION

1. Brush head
2. Water flosser nozzle
3. Hinged lid
4. Reservoir
5. Water hose
6. Water flosser pressure control dial
7. Water hose release button
8. Toothbrush ON/OFF button
9. Motion activated ring lights*
10. Water flosser ON/OFF button
11. Toothbrush speed control
12. Toothbrush speed indicator Lights
13. Toothbrush charge indicator lights
14. Power cord
15. Charging base
16. Wireless connection light
17. Magnetic brush head holder*

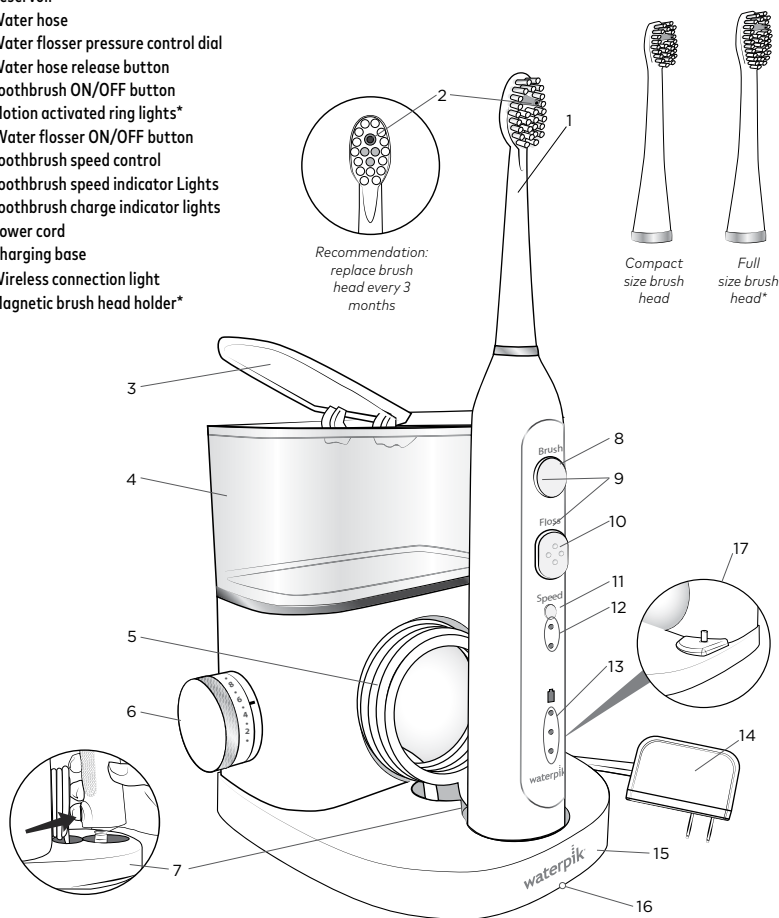Replacement brush heads may be purchased online at www.waterpik.com, or by phone from Water Pik, Inc. at 1-800-525-2774 (USA). To find detailed information please visit our website at www.waterpik.com.

*select models only

PREPARING THE UNIT / BRUSHING AND FLOSSING

PREPARING THE UNIT

Powering the Unit

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.

Attaching the Brush Head to the Power Handle

Position the brush head on the power handle as shown. Press firmly until it clicks into place. To remove, firmly pull the brush head to separate it from the power handle.

Attaching the Power Handle to the Hose

Position the power handle onto the hose in the base. Press the power handle firmly onto the hose attachment until it clicks into place. To remove, simply press the round button on the back side of the power handle.

Charging the Sonic-Fusion® 2.0 Power Handle

For a full charge, leave the handle on the charging base for at least 4 hours. If the toothbrush is to be stored for an extended period of time (six months or more), be sure to fully charge it (at least 4 hours) prior to storage.

Charge Lights

Toothbrush charge indicator lights will indicate level of battery charge for approximately 10 seconds after removal from the charging base and after use, and will then turn off. If only the bottom light is lit, then recharging is recommended to ensure uninterrupted use.

Wireless Connection Light








Blue light indicates the flosser base unit is communicating via Wireless Communication and is "Ready" for activation via FLOSS button.







Preparing the Reservoir

Remove reservoir and open hinged lid. Fill the reservoir with warm water. Close the lid and replace the reservoir on the base and press firmly.

Prime Unit (First Time Use Only)

Point tip into sink. Press FLOSS button to turn water on. Rotate pressure knob repeatedly between high and low until water flows. After water flows from the nozzle, press the FLOSS button to turn water OFF.

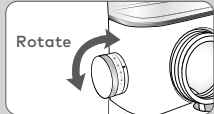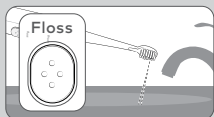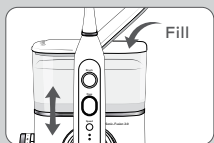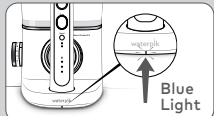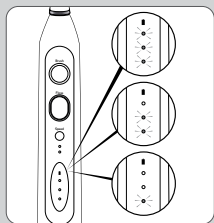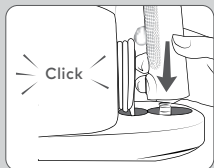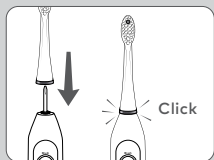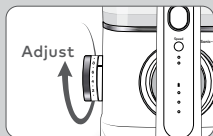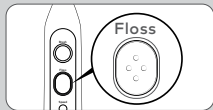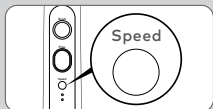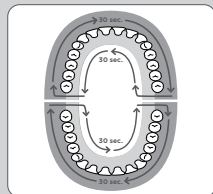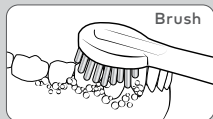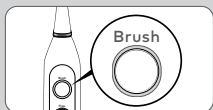BRUSH

Place the brush head in your mouth before pressing the **BRUSH** button.

Recommended Technique

For best results, position brush along gumline at a 45-degree angle. Press the **BRUSH** button to turn the brush ON. Use light pressure and move the toothbrush slowly allowing for the sonic action to clean each tooth individually. Brush can be easily disconnected from the hose if desired (See diagram on product description page).

Quad Timer

The Sonic-Fusion® 2.0 will create a short pause every 30 seconds to remind you to move to the next section of your mouth. You will brush the outside of the teeth for 30 seconds and then the inside for 30 seconds. This will be repeated on the other arch for a total of 2 minutes to clean all teeth.

2-Minute Timer

The brush will automatically shut off after 2 minutes of run time.

Brush Speeds

The Sonic-Fusion® 2.0 toothbrush gives you the option to choose between two brushing speeds: High and Low. If you wish to change to a different speed, you can do so by pressing the Speed button while brushing. The LED's below the Speed button indicate the current selection. High: 2 LEDs ON, Low: 1 LED ON.

FLOSS

Always place the brush head in your mouth prior to pressing the FLOSS button. Lean over the sink and slightly close lips enough to prevent splashing while still allowing water to flow from mouth into sink.

Adjusting the Pressure Setting

Turn the water flosser pressure control dial on the base of the unit to the lowest setting (pressure setting No. 1) for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

Recommended Technique


For best results, start with the back teeth and work toward the front teeth. Direct the brush head so that the water flosser nozzle is angled slightly toward the gumline. Move the brush head along the gumline pausing briefly between each tooth to allow water to clean between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth. When finished, empty any liquid left in the reservoir.

NOTE: If you used mouthwash or any other solution besides water, see "Using Mouthwash and Other Solutions" section.

BRUSH + FLOSS

To Brush + Floss at the same time, place brush head in mouth. Press **BRUSH** button, then press **FLOSS** button.


CARE AND MAINTENANCE

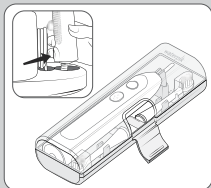Travel and Hose-Free Brushing

Toothbrush can be disconnected from the Sonic-Fusion® 2.0 for travel or hose-free brushing. Place the handle and desired brush head(s) in the provided case for protection when traveling.

NOTE: To disconnect the toothbrush from the water hose, press the round button on the back side of the power handle. See Product Description section.

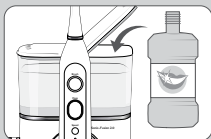Using Mouthwash and Other Solutions


Your Waterpik® Sonic-Fusion® 2.0 can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with brush head pointed into sink.

CARE AND MAINTENANCE

Before cleaning, unplug from the electrical outlet.

After Each Use

1. After each use, when all functions are turned OFF, remove the brush head from the handle and rinse brush head under running water.
2. Wipe the handle dry after each use.

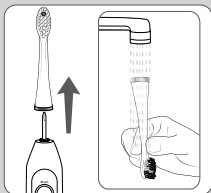General Cleaning / Removing Hard Water Deposits

- Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.
- Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 2 tablespoons of white vinegar to a full reservoir of warm water. Point the handle and brush head into sink. Turn Water Flosser ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water. This process should be done every 1 to 3 months to ensure optimal performance.

Leakage

If leaking due to a damaged hose, replacement hoses can be purchased at www.waterpik.com or contact Water Pik, Inc. customer service for a self service repair kit.

SONIC-FUSION TOOTHBRUSH RESET INSTRUCTIONS

In the event that your Sonic-Fusion 2.0 becomes ON/OFF unresponsive, there are two ways to reset your water flosser/toothbrush.

Radio Reset Steps:

1. Unplug the flosser base unit, wait 1 minute, Re-Plug base unit to AC-Mains. This performs a base radio reset.
2. Press and hold the Toothbrush BRUSH button. Then sequentially Press and hold the SPEED mode button; continue holding both buttons simultaneously. The toothbrush may run for 3 seconds and then automatically shut off at the 3 second interval.
3. Continue to HOLD BOTH buttons for 5 seconds.
4. After 6 seconds RELEASE BOTH BUTTONS. This shall reset the toothbrush to factory defaults. The toothbrush fuel gauge will also reset back to empty.
5. Place toothbrush in charging dock allowing radios to re-connect and begin a charging cycle. Note: the toothbrush will not run until adequate battery charge is provided back to the toothbrush indicated with at least a single solid (non-flashing) fuel gauge light.

Toothbrush Reset Steps:

1. Press and hold the toothbrush FLOSS button for a duration of 10 seconds.
2. After 10+ seconds, release the FLOSS button. At this point the toothbrush is reset to factory defaults. The toothbrush fuel gauge will also reset back to empty.
3. Place toothbrush in charging dock allowing radios to re-connect and begin a charge cycle. Note: the toothbrush will not run until adequate battery charge is provided back to the toothbrush indicated with at least a single solid (non-flashing) fuel gauge light.





Battery Removal

This product contains a rechargeable Lithium-Ion battery. **The battery is not replaceable.** Please do not dispose of the unit in the household waste. Remove the batteries per the instructions below and recycle at your local recycling center.

Note: This process is not reversible.

Remove the battery only at the end of the product's life for recycling. Doing so earlier will destroy the product. Deplete the battery charge until the brush no longer operates. Observe basic safety precautions to protect your eyes, hands, fingers, and work surface. Remove the battery as shown and recycle them and the rest of the product according to local regulations.

1. Run down the battery by operating the brush button until the motor does not run.
2. Obtain a towel or cloth, hammer, flat-head screwdriver, and wire cutters.
3. Remove and discard the brush head.
4. Cover the handle entirely with the towel or cloth.
5. Clamp the top of the handle against your work surface with one hand and strike the handle housing around the bottom perimeter on all sides to release the end cap latch. Continue striking and rotating the handle until the end cap releases.
6. Remove the internal components by turning the handle upside down and while holding on to the handle housing, press the shaft down firmly on a hard surface.
7. The internal components should slide out when the handle is positioned upside-down. If this does not occur, repeat steps 5 and 6 until the components release.
8. Remove the end cap from the internal assembly.
9. Slide the battery out of the internal assembly.
10. Cut the battery wires one at a time with the wire cutters.
11. Recycle and dispose of the battery properly.
12. The remaining product pieces should be discarded appropriately.

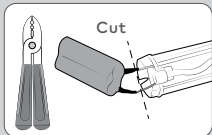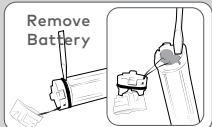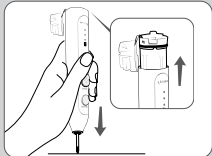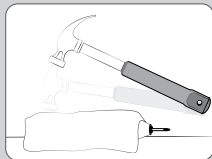LIMITED THREE YEAR WARRANTY

Service Maintenance

Waterpik® Sonic-Fusion® 2.0 has no consumer-serviceable electrical items and does not require routine service maintenance. For all accessories go to www.waterpik.com. You may also call our toll-free customer service number at 1-800-525-2774. Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Still have questions? Go to www.waterpik.com.

Water Pik, Inc. warrants to the original consumer of this product that it is free from defects in materials and workmanship for 3 years from the date of purchase.

Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace, at its discretion, the product which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. This limited warranty excludes accessories or consumable parts such as tips and brush heads. This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

For USA Only

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Limitation of Liability

This limited warranty is your sole remedy against Water Pik, Inc. and it shall not be liable for any damages, whether direct, indirect, incidental, special, consequential, exemplary or otherwise, including but not limited to property damage, arising out of any theory of recovery, including statutory, contract or tort.

Arbitration

BINDING ARBITRATION AND WAIVER OF RIGHT TO ASSERT MATTERS IN COURT AND RIGHT TO JURY TRIAL.

By purchasing, installing, or using this product you agree that any controversy or claim arising out of or relating to your purchase, use, modification, installation, or removal of this product shall be settled by binding arbitration. The arbitrator's decision will be final and binding. Judgment on the award rendered by the arbitrator may be entered in any court having jurisdiction over the matter. Further, you agree that you will not have the right to litigate any claim arising out of or in relation to your purchase, use, modification, installation, or removal of this product in any court or to have a jury trial on that claim.

Arbitration of any matters shall be administered by the American Arbitration Association and conducted by one arbitrator in accordance with the consumer arbitration rules of the American Arbitration Association in effect at the time of filing of the demand. The venue for all such arbitrations shall be Colorado and any such arbitration shall be subject to the Federal Arbitration Act and Colorado law.

WAIVER OF RIGHT TO PARTICIPATE IN CLASS ACTIONS OR PURSUE CLAIMS IN A REPRESENTATIVE CAPACITY.

Further, by purchasing, installing, or using this product you agree that arbitration must be on an individual basis. This means neither you, nor we, may join or consolidate claims in arbitration by or against other purchasers or users of Water Pik, Inc. products or the general public. Additionally, you may not litigate in court or arbitrate any claims as a representative or member of a class or in a purported representative capacity on behalf of the general public, other purchasers or users of this product, other persons or entities similarly situated, or in a private attorney general capacity.

For Canada Only

This device complies with Innovation, Science and Economic Development Canada (ISED)'s license-exempt RSSs. Operation is subject to the following two conditions: (1) This device may not cause interference; and (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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**Water Pik, Inc., a subsidiary of
Church & Dwight Co., Inc.**
Made in China. Designed and tested in the U.S.A.