

Register and Save!

Save over \$50!

Register at [Waterpik.com/Registration](https://www.waterpik.com/Registration) to receive exclusive discount code.

We'd Love to Hear from You!

Connect with Us
@Waterpik



Call: 1-800-525-2774
www.waterpik.com

Waterpik® 14 Day Challenge App

Make water flossing a habit! Try our new mobile App for the best water flossing experience.



Offer includes
50% off
your next order
of assorted
tips + tip
storage case.



waterpik® **FOR KIDS**
waterflosser

Congratulations!

You are now on your way
to better oral health!

✓ Register and Save

Register your product at [waterpik.com/registration](https://www.waterpik.com/registration) to receive over **\$50** in savings.

✓ Write a Review

Creating innovative Waterpik® products that inspire better oral health is our passion. We invite you to leave a review by sharing your thoughts on the Waterpik® Kids Water Flosser on [waterpik.com](https://www.waterpik.com) or the retailer's website.

✓ We're Here to Help

We hope you are loving your new Waterpik® Kids Water Flosser. If for any reason you are not satisfied or have any product questions, please contact our fast and friendly Customer Service Team based in Fort Collins, Colorado. Waterpik® Kids Water Flosser is backed by a 2-year manufacturer's warranty.

Call: 1-800-525-2774
Monday-Friday
7 am - 5 pm MST



Chat or E-mail
service@waterpik.com




**Quick Start
Guide Inside!**




GETTING STARTED

- 1 Fill reservoir with warm water and replace it firmly on the base.

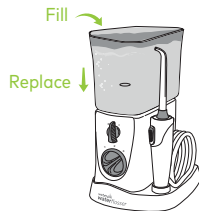


- 2 Insert tip firmly until it clicks into place. To remove tip, push the tip eject button.

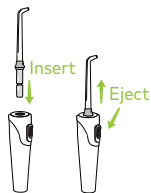

PRIMING UNIT (First Time Use Only)

- 3 Turn pressure to high, point tip into sink, turn power ON until water flows. Turn power OFF.

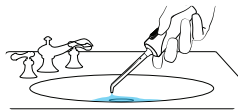

WATER FLOSS

- 4 Place tip in mouth. Start on low pressure. Turn power ON and increase pressure to your liking.




- 5 Lean over sink. Close lips enough to prevent splashing, while still allowing water to flow from mouth into sink.



TECHNIQUE

- 6 Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

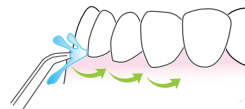

NEED HELP?

Visit [Waterpik.com/WelcomeWF](https://www.waterpik.com/WelcomeWF)

To watch a How To Video and find additional product support.

Call Us! 1-800-525-2774