

waterpik®
waterflosser™

Congratulations!

You are now on your way
to better oral health!

✓ Register and Save

Register your product at waterpik.com/registration to receive over \$50 in savings (includes 50% off your next order of assorted tips + tip storage case)

✓ Write a Review

Creating innovative Waterpik® products that inspire better oral health is our passion. We invite you to leave a review by sharing your thoughts on your new Water Flosser on waterpik.com or the retailer's website.

✓ We're Here to Help

We hope you are loving your new Water Flosser. If for any reason you are not satisfied or have any product questions, please contact our fast and friendly Customer Service Team based in Fort Collins, Colorado. your new Water Flosser is backed by a Two Year manufacturer's warranty.

Call: 1-800-525-2774 (USA)
Monday – Friday
7 am – 5 pm MST

Chat or E-mail
service@waterpik.com



**Quick Start
Guide Inside!**



Register and Save!

Save over \$50!

Register at Waterpik.com/Registration to receive exclusive discount code.

Offer includes
50% off
your next order
of assorted
tips + tip
storage case.



We'd Love to Hear from You!



Connect with Us
@Waterpik



Call: 1-800-525-2774 (USA)
www.waterpik.com

Waterpik® 14 Day Challenge App



Make water flossing a habit! Try our mobile App for the best water flossing experience.



Download on the
App Store




GET IT ON
Google play

Quick Start Guide

GETTING STARTED


- 1 With the unit turned **OFF(O)**, connect the charging cord to user supplied *USB charger. Plug the charger into the wall outlet, then connect the other end of the charging cord to the unit, and charge.

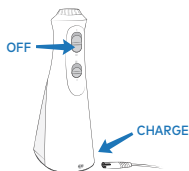

Charge Unit for 24 Hours Prior to First Use.

*USB charger requirements: USB Type-A, Rated Output: 5.0 VDC, 1.0 A


- 2 Press tip firmly into unit until it clicks into place. To remove tip, push the tip eject button.

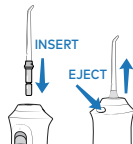

Caution:
Do not press tip eject button while unit is running.


- 3 Unplug charging cord from unit prior to filling with water and during use. Lift the flip top on reservoir and fill with lukewarm water.

Refill as needed.

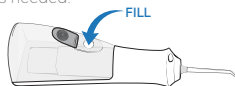

WATER FLOSS

- 4 For first time use, start on the **LOW(L)** setting.

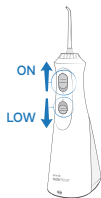

- 5 Lean over sink and close lips enough to prevent splashing while still allowing water to flow from mouth into the sink.



Turn unit **ON(O)**



Keep unit upright during use to ensure continuous water flow.

TECHNIQUE

- 6 Aim water at the gumline at a 90 degree angle.

Follow the gumline and pause briefly between teeth.



TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

NEED HELP?

Visit Waterpik.com/WelcomeWF
To watch a How To Video and find additional product support.

Call Us! 1-800-525-2774 (USA)