
waterpik®

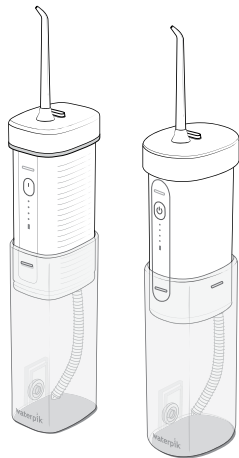Waterpik® Water Flosser Models WF-16 WF-17

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waterpik.com

SAFEGUARD STATEMENTS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not handle the user supplied USB wall charger with wet hands.
- Do not place or drop the user supplied USB wall charger into water or other liquid.
- Check the charging cord and user supplied USB wall charger for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. at 1-800-525-2774 (US) 1-888-226-3042 (Canada) or visit www.waterpik.com (US) www.waterpik.ca (Canada).
- Keep the charging cord and user supplied USB wall charger away from heated surfaces.
- Never drop or insert any object into any opening or hose.

- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not expose the battery or appliance to excessive temperatures.
- A user supplied USB-A wall charger (power adapter) with the following specifications is required to charge this device. **Rated Output : 5.0 VDC, 1.0 A.**
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Unplug USB wall charger when not in use.
- Do not use the Water Flosser for more than five minutes in each two hour period.
- Do not direct water from the Water Flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dental professional before using this product or any other oral hygiene aid.
- Fill reservoir with warm water or other dental professional recommended solutions only.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with a new Pik Pocket™ tip.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- There are no consumer serviceable electrical parts in this product and it does not require maintenance.
- Remove charging cord from the unit and replace charging port cover prior to filling the reservoir.
- If the supply cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.

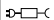
SAVE THESE INSTRUCTIONS

TABLE OF CONTENTS

Product Description	5
Getting Started	6
Cleaning and Troubleshooting Your Waterpik® Water Flosser	9
Battery Care, Disposal/Recycling, Battery Removal	10
Limited Two-Year Warranty	11


SYMBOLS

The following symbols may appear on your product and/or packaging. Some of these symbols may not be relevant in your region and are listed for information purposes only.

Symbol	Definition	Symbol	Definition	Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Detachable Power Supply		Do not use this device in a shower, bathtub and do not immerse in any liquid.		Regulatory Compliance Mark
	Manufacturer		Direct Current (DC) Voltage		Authorized Representative in the European Community		Separate collection for waste electrical and electronic equipment.
IPX7 or IPX4	Ingress Protection Code		Alternating Current (AC) Voltage		CE marking of conformity		
	IEC Class II Equipment		Do not immerse this device in any liquid.		UK Conformity Assessed		

PRODUCT DESCRIPTION

1. Power/Mode button
2. LED mode indicators
3. Battery charge indicator
4. Tip
5. Tip eject button
6. Reservoir
7. Reservoir fill door
8. Charging Port Cover

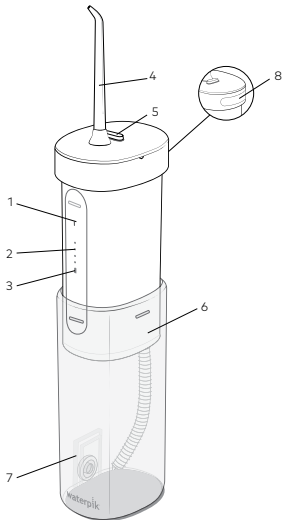TIPS

- A. Orthodontic Tip***
- Braces/General use
- B. Pik Pocket™ Tip***
- Periodontal Pockets/Furcations
- C. Plaque Seeker™ Tip***
- Implants/Crowns/Bridges/
Retainers/General use

D. Classic Jet Tip**

- General use


E. Tongue Cleaner**

- Fresher breath

F. Travel Plug

- Keeps excess water from
draining out while traveling

NOT ALL TIPS INCLUDED IN ALL MODELS



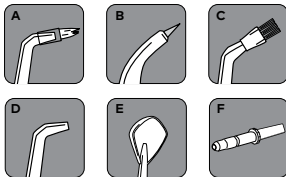
*Replace Every 3 months **Replace Every 6 months

Replacement tips/attachments may be purchased online at www.waterpik.com (US) www.waterpik.ca (Canada), or by phone from Water Pik, Inc. at 1-800-525-2774 (US) or 1-888-226-3042 (Canada).

USB Charging Cable
20030405-2 White

GETTING STARTED

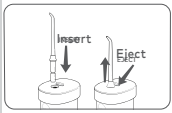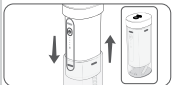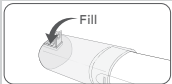Charging Unit Prior to First Use

Charge unit for approximately 4 hours prior to first use. Plug the user supplied USB wall charger into a wall outlet and connect the USB charging cable to the unit (see images). Refer to the *IMPORTANT SAFEGUARDS* section for required USB wall charger specifications.

Ensure contacts on BOTH the cord and the unit are clean and dry, and the reservoir is empty and collapsed. The LED light will emit a slow glow on/off when unit is charging. Once the unit is fully charged and the charger is connected the LED will remain on until the charger is disconnected. A full charge will last up to 4 weeks of typical use. When the batteries are in need of charging (approximately 30% charge left) the charge indicator light emit a rapid triple blink during and after use to indicate that you should recharge the unit.

A complete recharge will take 3-4 hours. If you use your Water Flosser once a day or less, it is recommended you charge every 3-4 weeks. If you use it twice a day or more, charge more frequently.

Put the protective charging port cover back into place before use. The charging port cover should only be removed for charging.

Filling the Reservoir


Remove the charging cable from the unit and replace charging port cover prior to filling reservoir and product use. Pull down on reservoir to open. Unit will "click" into place when fully extended. Lift flip top on reservoir and fill the reservoir with lukewarm water.

Collapse by pushing reservoir up. Reservoir will "click" into place when fully closed.

Inserting and Removing Tips

Insert tip at the top of the unit. Press firmly until the tip clicks into place. To remove the tip from the unit, press the tip eject button and pull tip from the unit.

CAUTION: Do not press tip eject button while unit is running.

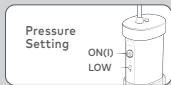Adjusting the Pressure Setting

1. 3 modes – low, med, high pressure settings.
2. Unit will remember last setting when turned off (press and hold if within 5s of turn on, short press after running for more than 5s).
3. Unit will cycle through the modes when the button is repeatedly pressed within 5s of turn on.


Recommended Technique

Lean low over sink and place tip in mouth. Aim the tip toward gumline, and turn the unit ON.


For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

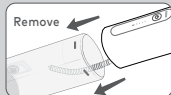
When Finished

Turn the unit OFF. Empty any liquid left in the reservoir. If desired, the reservoir may be removed from the power handle by simply pulling down on the reservoir to remove.

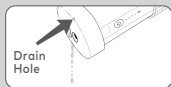Drain unit after use to prevent residual water from leaking out. Tilt unit on its side and allow water to drain from the drain holes at the top.

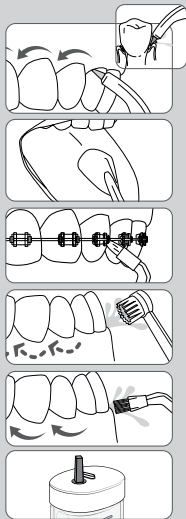NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket™ tip, set the unit to the lowest pressure setting. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.

Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure.

Increase pressure as you prefer.

Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

Toothbrush Tip

Place the toothbrush tip in mouth with the brush head on the tooth at the gumline. The toothbrush tip can be used with or without toothpaste. Using a light pressure (bristles should not bend) massage the brush back and forth with very short strokes – much like you would with a manual toothbrush.

Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.


Travel Plug

When travelling, place travel plug in unit to keep excess water from draining out.

NOTE: Do NOT run unit while travel plug is in unit.

Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. For ease of adding other solutions remove the reservoir (see "When Finished" section) and add from opening on reservoir. If you are adding additional water you may replace the reservoir and fill reservoir fully with the fill door. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

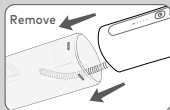Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild nonabrasive cleanser to wipe exterior. The water reservoir is removable for easy cleaning. It is recommended that the reservoir is removed and cleaned weekly.

Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water. This process should be done every 1 to 3 months to ensure optimal performance.

BATTERY CARE, DISPOSAL/RECYCLING, BATTERY REMOVAL

Battery Care

If the Water Flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage. Do not expose the battery or appliance to excessive temperatures.

Disposal/Recycling

This product contains a **rechargeable lithium ion battery**. The battery is not replaceable. At the end of product life, do not dispose of in household garbage or recycling bins. Please recycle at your nearest recycling center or take to hazardous waste collection point to comply with local and state regulations.

Please contact your local waste administrative department for any questions.


Note: The rechargeable battery is not replaceable. Remove only at end of product life. This process is not reversible and voids the product warranty.

Battery Removal

To remove the rechargeable battery for recycling, please discharge the battery fully before removing.

Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and surface on which you work.

1. The power should be off and the battery fully discharged before removal. It may take several power cycles before the battery is fully discharged.
2. Cut off top of unit.
3. Slide battery out of chassis. Cut wires individually (avoid cutting positive and negative at the same time to prevent shorting).

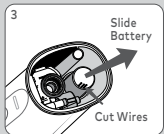Service Maintenance

Waterpik® Water Flossers have no consumer serviceable electrical items and do not require routine service maintenance. For all parts and accessories go to www.waterpik.com (US) or www.waterpik.ca (Canada). You may also call our toll-free customer service number at 1-800-525-2774 (US) or 1-888-226-3042 (Canada).

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Still have questions? Go to www.waterpik.com (US) or www.waterpik.ca (Canada).

Limited Two Year Warranty

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. This limited warranty excludes accessories or consumable parts such as tips, etc. Use with improper voltage will damage product and VOID warranty.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

FOR USA ONLY

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by Water Pik, Inc. could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Limitation of Liability

This limited warranty is your sole remedy against Water Pik, Inc. and it shall not be liable for any damages, whether direct, indirect, incidental, special, consequential, exemplary or otherwise, including but not limited to property damage, arising out of any theory of recovery, including statutory, contract or tort.

Arbitration

BINDING ARBITRATION AND WAIVER OF RIGHT TO ASSERT MATTERS IN COURT AND RIGHT TO JURY TRIAL.

By purchasing, installing, or using this product you agree that any controversy or claim arising out of or relating to your purchase, use, modification, installation, or removal

of this product shall be settled by binding arbitration. The arbitrator's decision will be final and binding. Judgment on the award rendered by the arbitrator may be entered in any court having jurisdiction over the matter. Further, you agree that you will not have the right to litigate any claim arising out of or in relation to your purchase, use, modification, installation, or removal of this product in any court or to have a jury trial on that claim.

Arbitration of any matters shall be administered by the American Arbitration Association and conducted by one arbitrator in accordance with the consumer arbitration rules of the American Arbitration Association in effect at the time of filing of the demand. The venue for all such arbitrations shall be Colorado and any such arbitration shall be subject to the Federal Arbitration Act and Colorado law.

WAIVER OF RIGHT TO PARTICIPATE IN CLASS ACTIONS OR PURSUE CLAIMS IN A REPRESENTATIVE CAPACITY.

Further, by purchasing, installing, or using this product you agree that arbitration must be on an individual basis. This means neither you, nor we, may join or consolidate claims in arbitration by or against other purchasers or users of Water Pik, Inc. products or the general public.

Additionally, you may not litigate in court or arbitrate any claims as a representative or member of a class or in a purported representative capacity on behalf of the general public, other purchasers or users of this product, other persons or entities similarly situated, or in a private attorney general capacity.

FOR CANADA ONLY

This device complies with Industry Canada's license-exempt RSSs. Operation is subject to the following two conditions: (1) This device may not cause interference; and (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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©2023 Water Pik, Inc.
Church & Dwight Co., Inc.**
Form No. 20032326-F AB